

IMPORTANT SAFETY INSTRUCTIONS

When using an electric massager, basic precautions should always be followed, including the following:

Read all instructions before using this massager.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) Do not use while bathing or in a shower.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) This massager is not intended for use by, on, or near children or people with certain disabilities.
- 4) Use this massager only for its intended use as described in this manual.
- 5) Keep the massager away from heated surfaces.
- 6) Never drop or insert any object into the massager.
- 7) To disconnect, turn all controls to the off position, then remove batteries from massager.
- 8) Do not use in bed.
- 9) Never use while sleeping.
- 10) Never operate this massager if it is damaged, if it is not working properly, or if it has been dropped into water. Return the massager to an authorized service center for examination and repair.
- 11) Do not use near any sharp objects to prevent massager from being punctured.
- 12) **CAUTION: CHOKING HAZARD.** This is not a toy. Not to be used as a baby or infant pillow. Do not open this massager for any reason. If this massager is damaged or opened, discard immediately out of reach of children.

SAVE THESE INSTRUCTIONS OPERATING INSTRUCTIONS

Insert Batteries

Before using your massager you must insert 2 AA batteries (not included). Unzip cover to access the battery compartment. Insert batteries. Align positive (+) and negative (-) as indicated inside the battery compartment. Close zipper.

Note: Use only alkaline batteries. DO NOT MIX old and new batteries. DO NOT mix alkaline, standard (carbon-zinc), or rechargeable (nickel cadmium) batteries. Remove batteries when storing massager for an extended period of time.

Your Portable, Battery Operated Massaging Seat Cushion is designed to be used in your home or office. Elastic straps are attached to the seat, which allow you to securely fasten the massager to almost any chair.


This massager has been designed for ease of operation. There are 3-push button switches on the control panel to control the vibrating motors:

1st switch	UPPER BACK
2nd switch	LOWER BACK LEFT
3rd switch	LOWER BACK RIGHT

Each switch has an OFF position. Try each motor setting to create the custom message you desire.

INSTRUCTIONS FOR USING BATTERIES

- Always purchase the correct size and grade of battery most suitable for the intended use.
- Replace all batteries of a set at the same time.
- Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (nickel cadmium) batteries.
- Clean the battery contacts and also those of the device prior to battery installation.
- Ensure the batteries are installed correctly with regard to polarity (+ and -).
- Remove batteries from equipment which is not to be used for an extended period of time.
- Remove used batteries promptly.

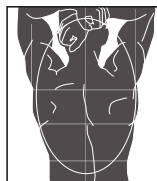

DISPOSAL OF BATTERIES

When replacing the batteries, make sure you dispose of the old batteries in accordance with applicable regulations. They should be disposed of in a location specifically designated for that purpose so they can be recycled safely and without posing any risk for the environment. Do not burn or bury them.

CAUTIONS

- Choking hazard. Keep out of reach of children. This is not a toy. Not to be used as a baby or infant pillow. Read instruction booklet before using.
- This massager should not be used in the following cases: Tuberculosis, benign and malignant tumors, hemorrhages, inflammations of the skin, phlebitis and thrombosis, open and fresh wounds, bruises, broken skin or varicose veins, undetermined calf pain, or joint pain.
- If you are undergoing treatment for cancer, check with your physician before using.
- **Consult your physician should you have any questions regarding use of this massager.**
- Do not use on any unexplained pain or swollen muscles, or following a serious muscular injury, before consulting your physician.
- Keep away from children.
- Discontinue use of massager if pain increases or swelling occurs after application.
- If you are pregnant, consult your physician before using.
- If you are taking any prescription medications, consult your physician before using. Individuals with certain disabilities should consult their physician before using.

MAINTENANCE INSTRUCTIONS

- 1) Store massager in a secure place away from excessive heat or sharp objects.
- 2) To clean the unit, remove the batteries and then wipe down with a damp cloth. Do not use any strong chemical cleaners such as benzene.

RECOMMENDED MASSAGING TIMES

Generally, massage should last 15 minutes on any one particular part of the body. Treatment may be applied once or twice daily as desired.

As metabolism varies among individuals, so does their response to massage. You will soon adjust the time to suit your particular needs and personal comfort. If you have any doubts regarding treatment, we recommend you consult your physician.

MASSAGING BACK & SEAT CUSHION

3 powerful motors provide soothing, vibrating massage to two massage zones

battery operated cord-free convenience (2 x AA not included)

easy-to-use push-button controls



upper back

lower back

soft-cushioned design provides extra comfort and flexibility

strapping system allows you to fasten to almost any chair

portable for home or office

LIMITED ONE YEAR WARRANTY

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and for compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure. In New Zealand our goods come with a guarantee that cannot be excluded under the Consumer Guarantee Act 1993.

Conair Australia Pty Ltd. will repair or replace (at our option) your appliance free of charge for twelve months from the date of purchase if the appliance is defective in workmanship or materials.

This warranty is an additional benefit and does not affect your legal rights. To obtain service under this warranty, please read the Warranty Leaflet included or contact Conair Conair Customer Service at aus_info@conair.com. To make a warranty claim you must retain your proof of purchase.

CONAIR®

©2018 Conair Australia Pty. Ltd.
Suite 101, 18 Rodborough Rd
Frenchs Forest NSW, 2086
Customer Service Hotline
Australia: 1800 650 263
New Zealand: 0800 266 247
E-mail: aus_info@conair.com
www.conairaustralia.com.au

IB-18/344

MASSAGING BACK & SEAT CUSHION

Portable, Battery Operated, provide soothing, vibrating massage



body benefits®
by CONAIR®

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Model CBMIRBA